

Randy Cole
Regenerative Research
Denver, CO 80210



Project 135319
Lab # 128699
Date Rec'd 9/22/2025
Report Issued 9/29/2025

Project

Sample PT-141 10 mg PT25-01

Certificate of Analysis

<u>Analyte</u>	<u>Result</u>	<u>LOQ</u>	<u>Units</u>	<u>Method</u>	<u>Date</u>	<u>% of Label</u>
Microbiology						
Endotoxins	ND	1.0	EU/mL	USP 85	9/29/2025	
Peptide Analysis						
Chromatographic purity	99.85	0.5	%	HPLC-UV/MS	9/23/2025	
Total peptide mass	11.0	0.5	mg	HPLC-UV/MS	9/23/2025	110.01
Mass spectral identification	Confirmed			HPLC-UV/MS	9/23/2025	

The data presented are from the analysis of the sample shown and meet Krause Analytical internal quality assurance criteria unless otherwise flagged.
Methods shown reference current Krause Analytical SOPs
ND - Not detected LOQ - limit of quantification
Copyright © 2025 Krause Analytical LLC. All rights reserved.

Respectfully submitted,


Mark C. Krause
Laboratory Director




Lab #: 128699 PT-141 10 mg PT25-01

Chromatogram

Response_

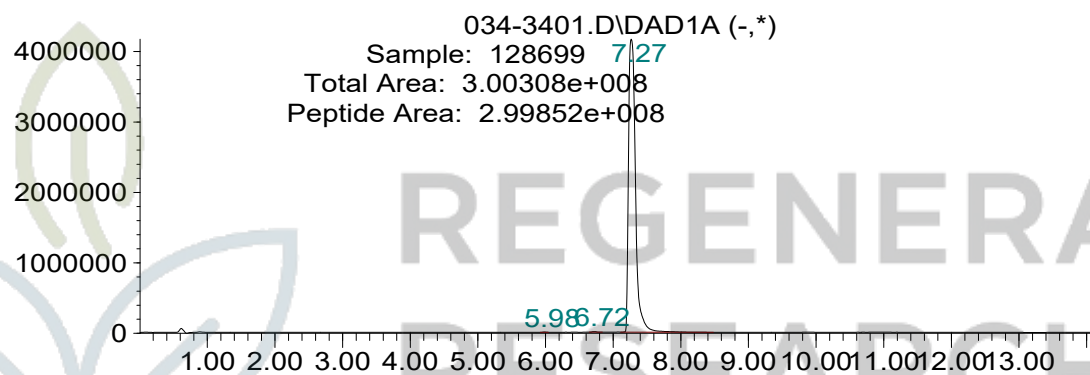


Response_

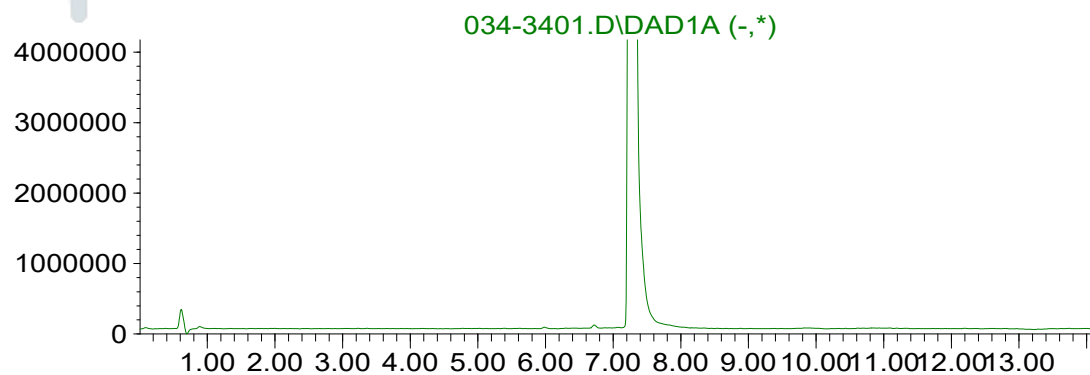

Time



Mass spectrum/Reference Spectrum

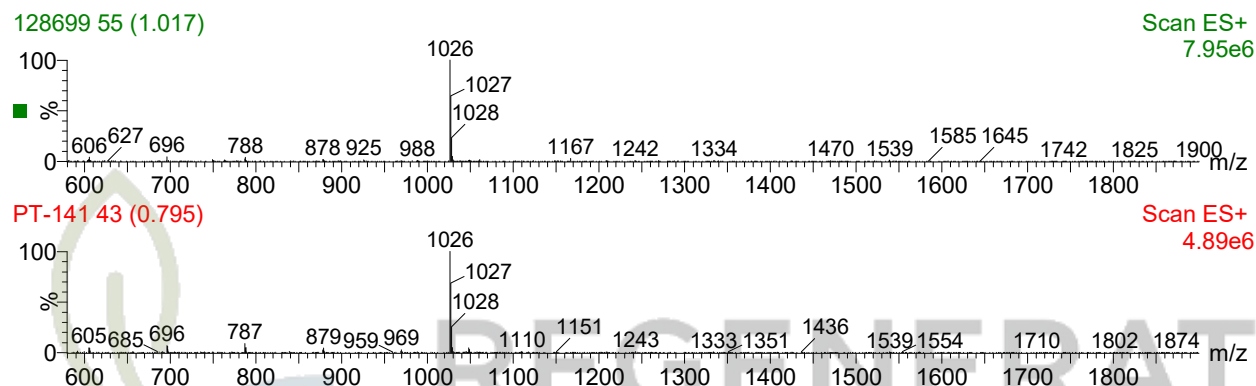

Method Summary

- 2 mL of purified water is added to the lyophilized powder in the vial, and the contents mixed to dissolve the lyophilized powder.
- An aliquot is taken from the vial and diluted to contain approximately 500 mg/L of the peptide.
- The diluted sample is analyzed by HPLC-UV-MS.
- The mass spectrum obtained is compared to an authentic standard of the peptide for identification.
- The total area of all of the peaks in the chromatogram is calculated, and the area of the peak of the peptide is divided by the total area to obtain the chromatographic purity value, reported in percent.
- The area of the peptide is compared to the area of the peptide peak in the known standard to obtain a concentration in the solution. This concentration is used to calculate the total mass of peptide in the vial, which is compared to the stated mass (label claim) and reported as both total mass in the vial and as a percent of the label claim.